

SELUTION SLR™ DEB

DRUG ELUTION DONE BETTER

SUBGROUP ANALYSIS

SELUTION DENOVO COMPLEX PCI¹

A pre-specified subgroup analysis of a prospective, multi-center, randomised, open-label, non-inferiority trial comparing a SELUTION SLR DEB strategy vs. systematic DES strategy in patients with complex PCI

COMPLEX PCI SUBGROUP 835 PATIENTS

**DEB
STRATEGY**
421 PATIENTS

1:1

RANDOMIZATION


**DES
STRATEGY**
414 PATIENTS

ASSESSING TVF AT 12 MONTHS

SUBGROUP ANALYSIS RESULTS

IN COMPLEX PCI, A **SELUTION SLR DEB STRATEGY**
DEMONSTRATES THE BENEFIT OF MINIMIZING STENTS

TVF AT 12 MONTHS*

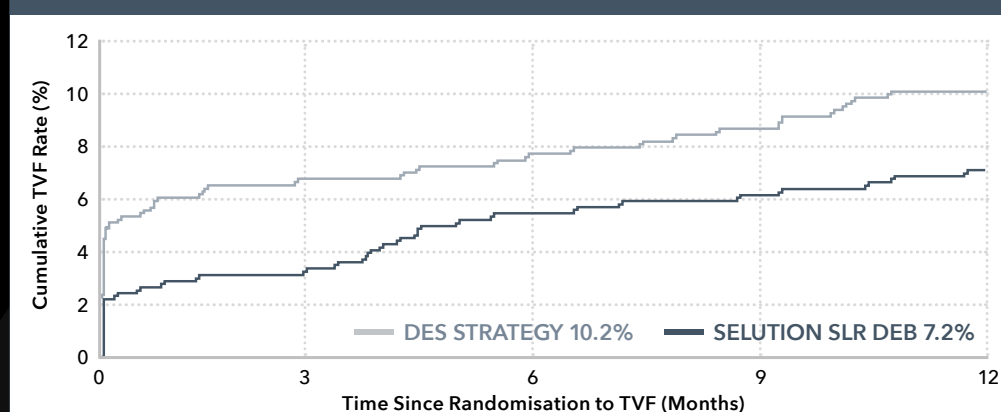

10.2%

TARGET VESSEL FAILURE
DES STRATEGY

VS

7.2%

TARGET VESSEL FAILURE
SELUTION SLR DEB

* Intention to Treat

COMPOSITE ENDPOINT OF TVF

63% LESS
CARDIAC DEATH

0.7% SELUTION SLR DEB VS. 1.9% DES STRATEGY

65% LESS
TV-MI

2.6% SELUTION SLR DEB VS. 7.5% DES STRATEGY


COMPARABLE
CD-TVR

4.3% SELUTION SLR DEB VS. 4.1% DES STRATEGY


SELUTION SLR DEB

SHOWED LOW ACUTE AND LATE SAFETY EVENT RATES

SAFETY ENDPOINTS IN COMPLEX PCI SUBGROUPS AT 1 YEAR

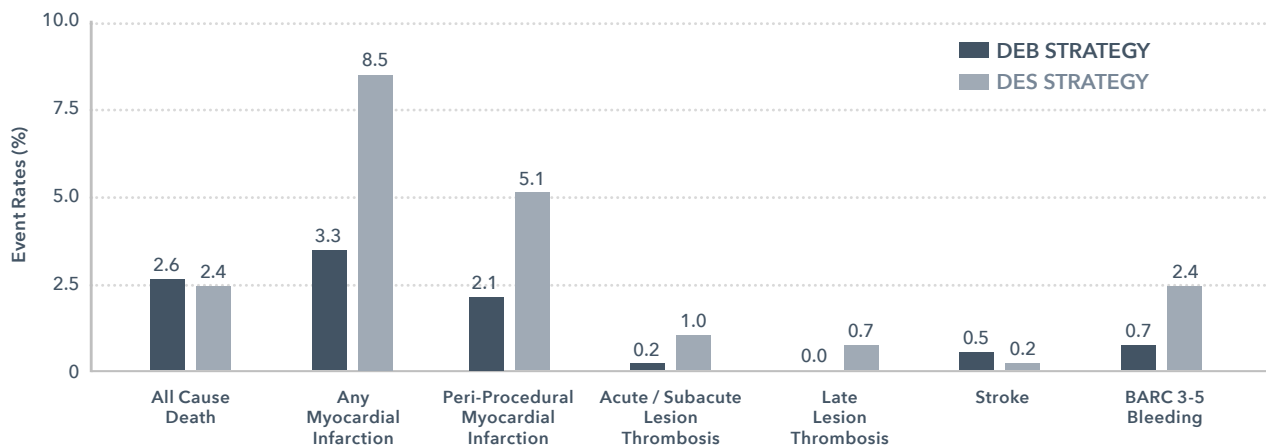


COMPLEX PCI CRITERIA AND CHARACTERISTICS

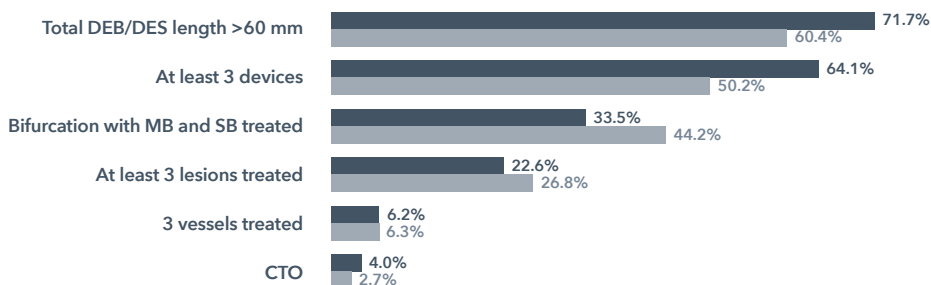

SELUTION SLR DEB

PROVIDES ENHANCED DELIVERABILITY WITH LESS FORCE REQUIRED COMPARED TO LEADING DCB COMPETITORS²

REQUIRED FORCE TO CROSS (MAX FORCE (G))

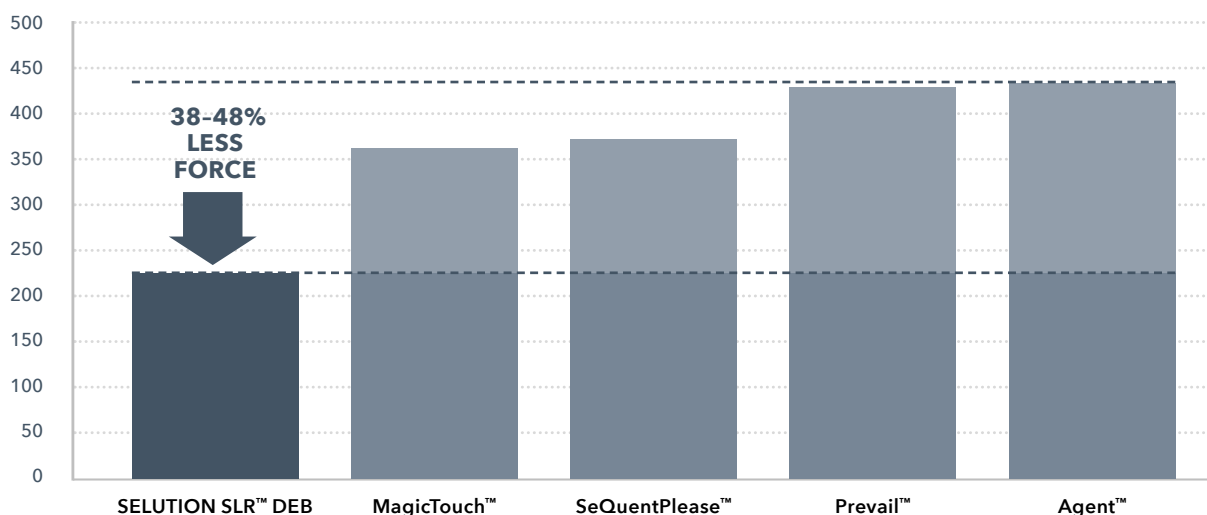

1. SELUTION DeNovo Trial Complex PCI Subgroup Analysis, Oral Presentation, EuroPCR 2026 This pre-specified subgroup analysis was not powered for statistical significance; observed differences are descriptive in nature.

2. Cordis 2026 Data on file.

For Healthcare Professionals Only.

Important information: Prior to use, refer to the instructions for use supplied with this device or available at www.cordislabeling.com for indications, contraindications, suggested procedure, warnings, and precautions. As part of its continuous product development policy, Cordis reserves the right to change product specifications without prior notification. Please contact your Cordis representative for additional product availability information. SELUTION SLR Drug-Eluting Balloon is manufactured by MedAlliance, LLC and its affiliates. SELUTION SLR is a trademark of M.A. Med Alliance SA, MedAlliance, LLC and M.A. Med Alliance SA are Cordis companies. The SELUTION SLR PTCA DEB fulfills the requirements for CE marking certified at Notified Body 0344. CORDIS, Cordis LOGO, are trademarks of Cordis and may be registered in the US and/or in other countries. All other marks are the property of their respective owners.

©2026 Cordis. All Rights Reserved. 100690689-2



Scan the QR Code to learn more.

