


*Cordis*<sup>™</sup>

# ADROIT<sup>™</sup>


## Guiding Catheter

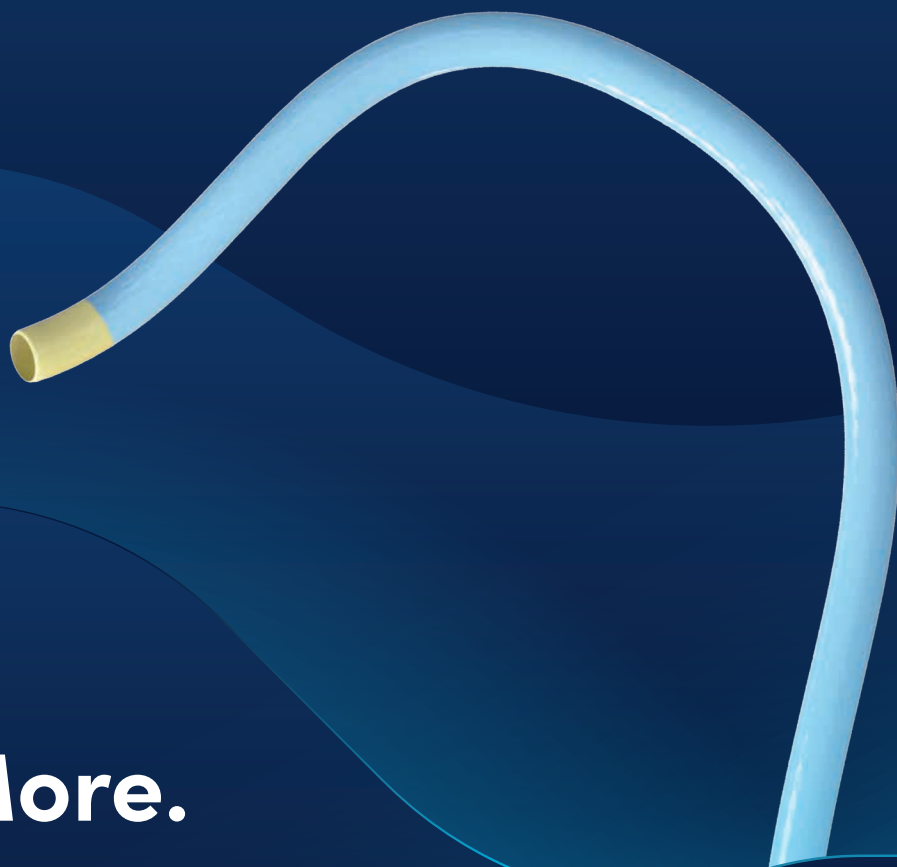

**Get More.**

Control  
Back-up support  
Space to work



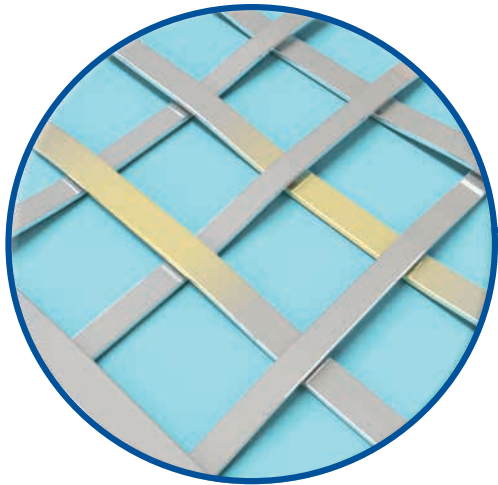


Hybrid design features two different sizes of wire braid

## BACK-UP SUPPORT

- Engineered with Hybrid braided wire technology
- Controlled Back-up Support

Back-up Support<sup>1,2</sup>

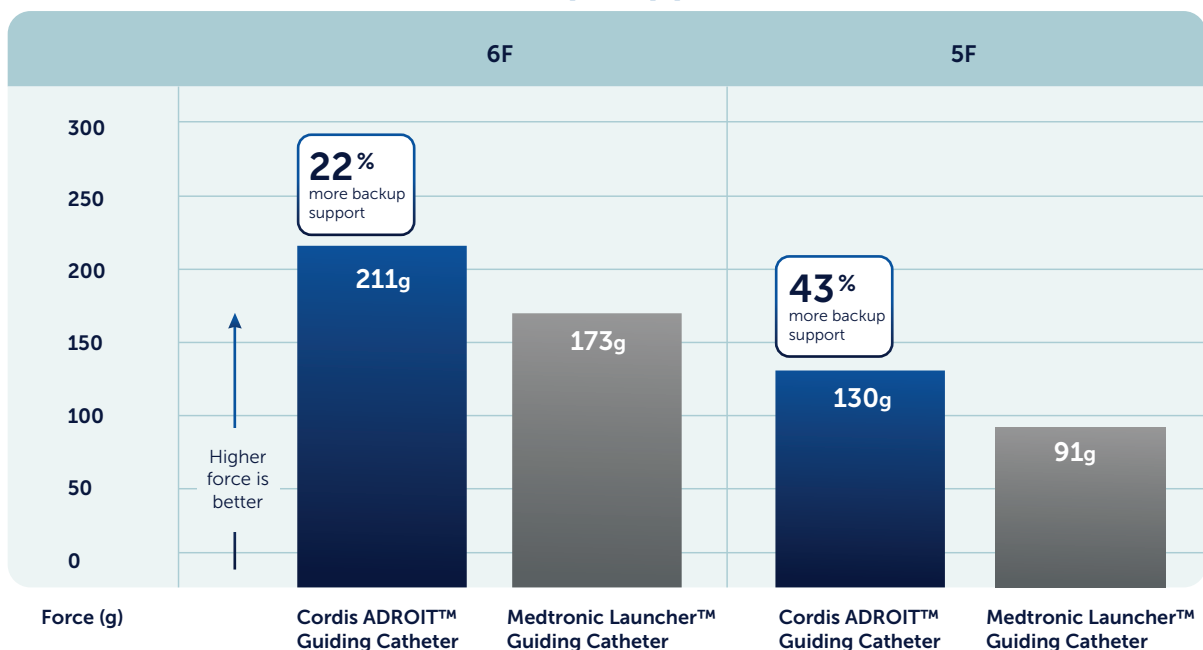

<sup>1</sup> Cordis 2012 and 2015 Data on file.


<sup>2</sup> For 6F configurations, sample size was  $\leq 3$ .

# CONTROL


## Gain greater positioning confidence through:

- Hybrid braided wire technology that delivers excellent torque response
- Multi-segment design that provides kink resistance and minimizes whip

### Torque Response<sup>1,2</sup>

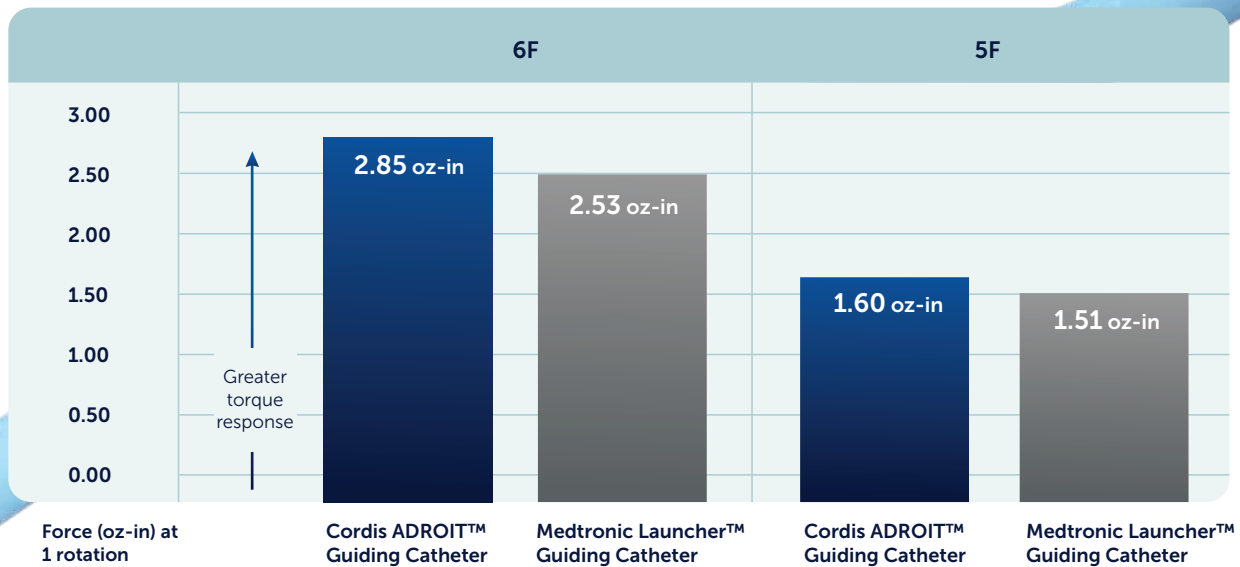

### Whip<sup>1,2</sup>

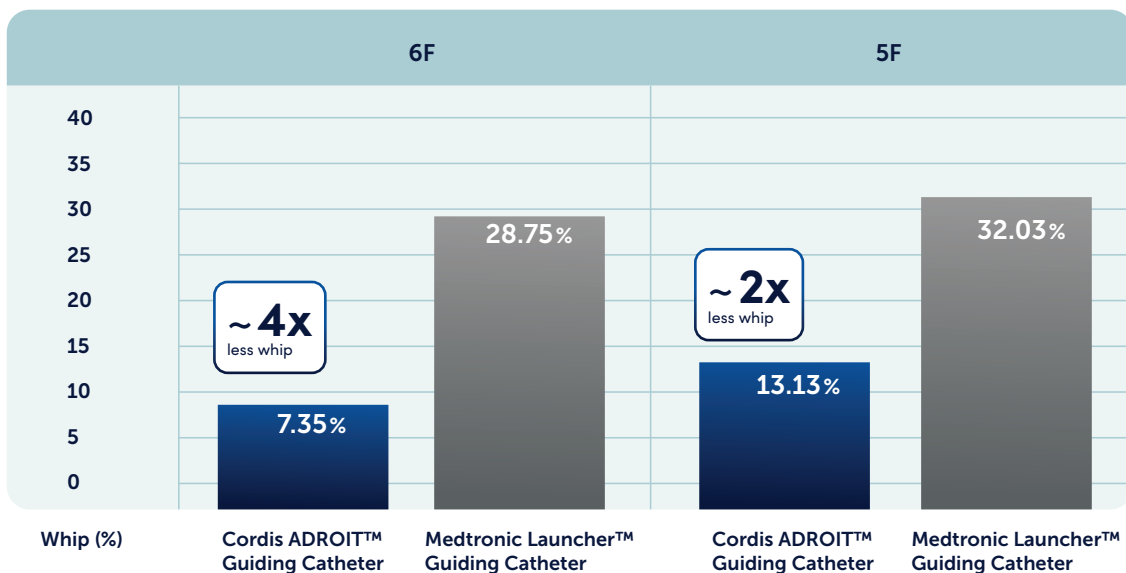


<sup>1</sup> Cordis 2012 and 2015 Data on file.

<sup>2</sup> For 6F configurations, sample size was ≤ 3.


# SPACE TO WORK

Experience more room to maneuver with the available Larger Lumen ID (6F)\*:

- .072" ID provides device compatibility and provides good visualization
- Large lumen to perform kissing balloon technique
- PTFE lining provides lubricious lumen for smooth delivery
- True lumen Technology ensures consistent ID from hub to tip

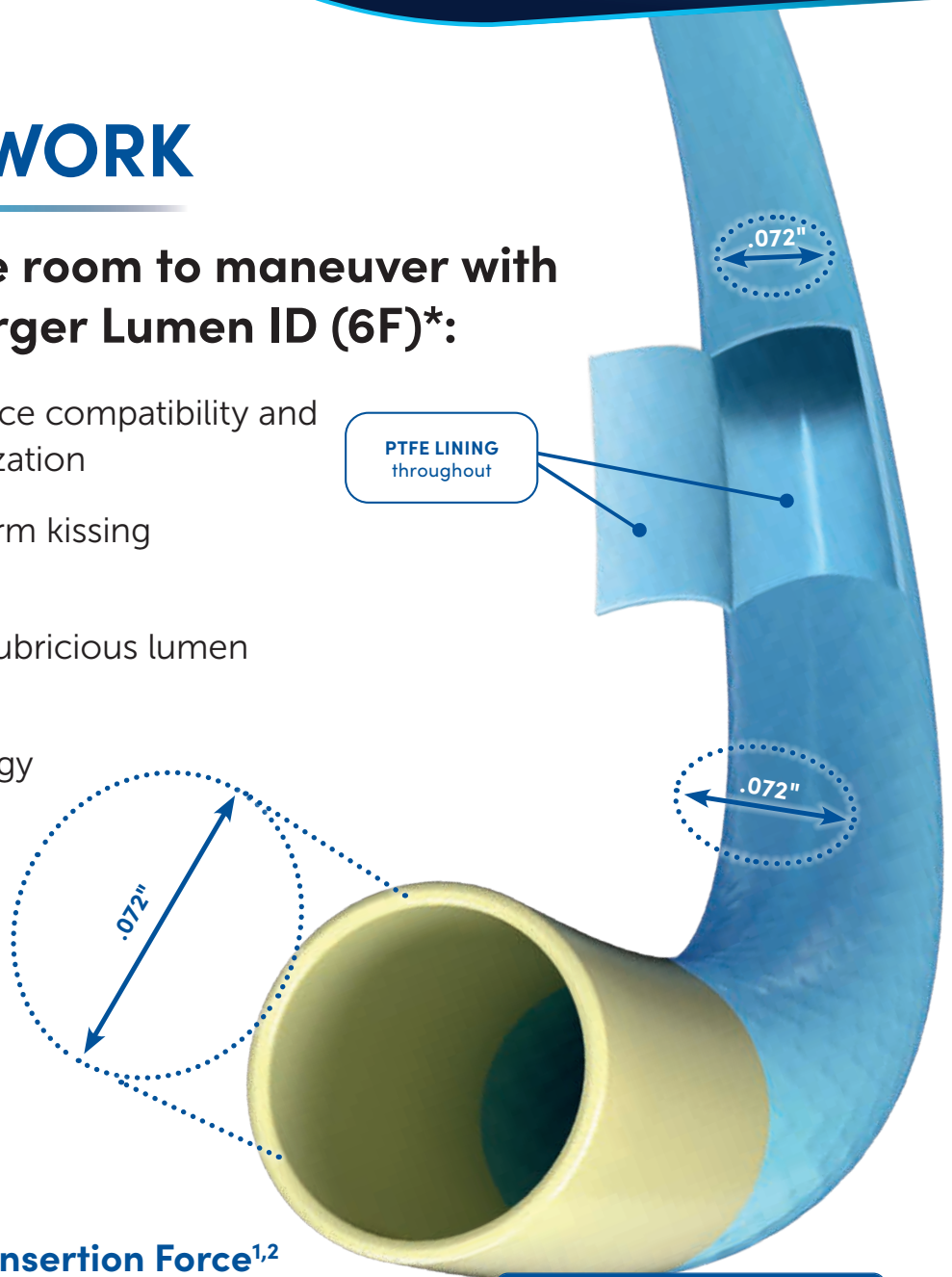

## Kissing Balloon Insertion Force<sup>1,2</sup>

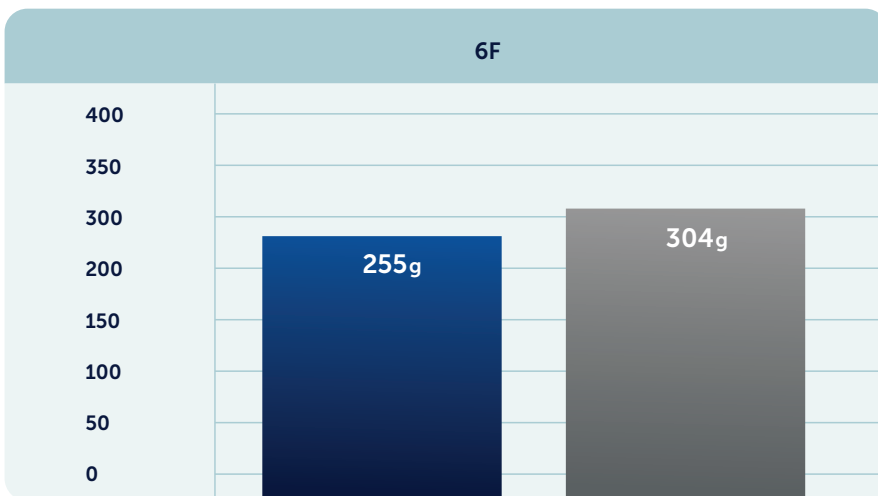


Force to track second balloon (g)

Cordis ADROIT™ Guiding Catheter

Medtronic Launcher™ Guiding Catheter

## Inner Lumen Diameter and Smoothness

Cordis ADROIT™ Guiding Catheter

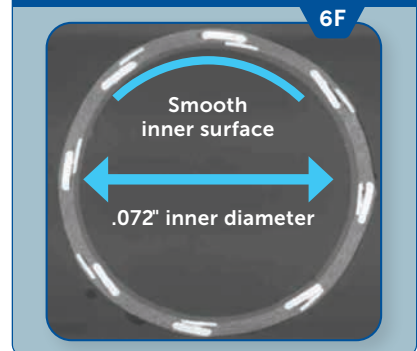

Hybrid braided wire technology enables large lumen diameter (.072")

\* Compared to Cordis VISTA BRITE TIP™ Guiding Catheter.


<sup>1</sup> Cordis 2012 and 2015 Data on file.

<sup>2</sup> For 6F configurations, sample size was ≤ 3.

# DESIGNED FOR ATRAUMATIC PLACEMENT

## Engineered to avoid vessel wall damage via:

- Longer, soft radiopaque tip provides enhanced visibility\* and allows active cannulation.
- A multi-segment design that provides kink resistance and minimizes whip

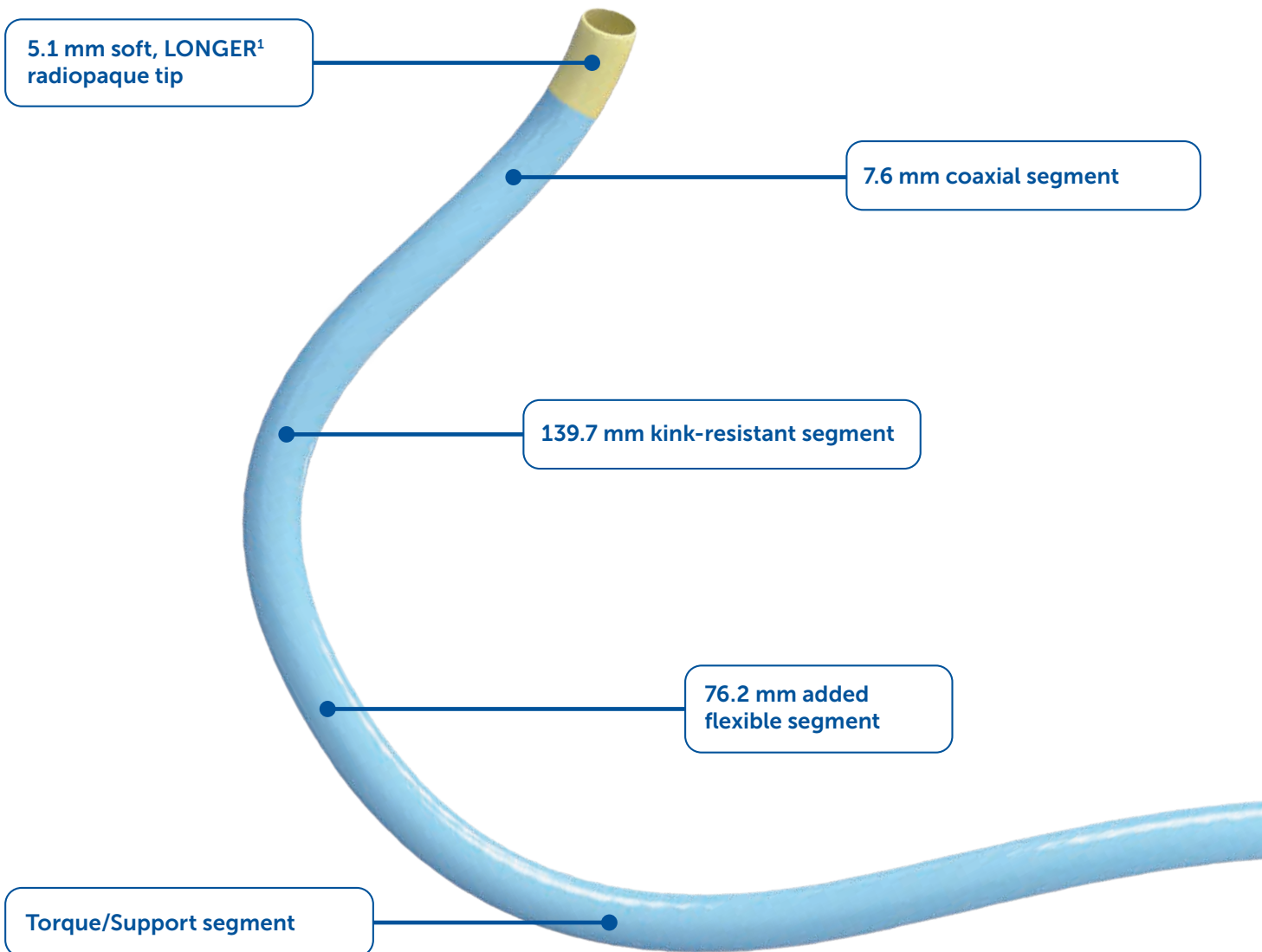

\* Compared to the Cordis VISTA BRITE TIP™ Guiding Catheter.

<sup>1</sup> The primary differences between the ADROIT™ Guiding Catheter and the VISTA BRITE TIP™ Guiding Catheter are the larger inner diameter of the ADROIT™ Guiding Catheter and the three transition segments in the ADROIT™ Guiding Catheter as opposed to two in the VISTA BRITE TIP™ Guiding Catheter.



ADROIT™ Guiding Catheter

# Get More.

## Control

A multi-segment design that provides kink resistance and minimizes whip

## Back-up support

Hybrid braid construction delivers support, torque control and lumen integrity


## Space to work

Large lumen to perform kissing balloon technique and PTFE lining provides lubricious lumen for smooth delivery

**Cordis**™

**Cordis ADROIT™** Guiding Catheter


A large portfolio of shapes and configurations are available to meet varying procedural and anatomical demands.

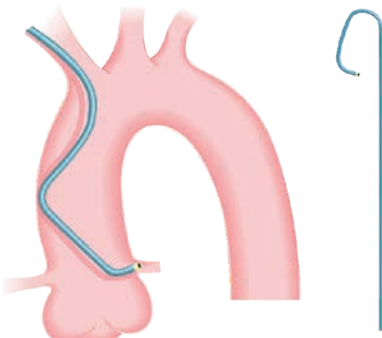
Left Coronary Artery	Shape	Description	PRODUCT CODE	
			ADROIT™ 6F	ADROIT™ 5F
 <p>Judkins Left</p>	JL 3		672-000-00	558-000-00
	JL 3.5		672-002-00	558-002-00
	JL 3.5 SH	2 Side Holes	672-003-00	
	JL 3.5 ST	Short Tip	672-012-00	
	JL 4		672-004-00	558-004-00
	JL 4 (90 cm)		672-004-90	
	JL 4 (125 cm)		672-040-25	
	JL 4 SH	2 Side Holes	672-005-00	
	JL 4 ST	Short Tip	672-014-00	
	JL 4 LBT	Long BRITE TIP™	672-004-0L	
	JL 4	ECOPAK (5 units)	672-004-0E	
	JL 4.5		672-006-00	558-006-00
	JL 4.5 ST	Short Tip	672-016-00	
	JL 5		672-008-00	558-008-00
JL 6		672-010-00		

 <p>Judkins Curved Left</p>	Shape	Description	ADROIT™ 6F	ADROIT™ 5F	
	JCL 3.5			672-026-00	558-026-00
	JCL 3.5 SH	2 Side Holes	672-027-00		
	JCL 4		672-028-00		

 <p>Amplatz Left</p>	Shape	Description	ADROIT™ 6F	ADROIT™ 5F	
	AL .75			672-034-00	558-034-00
	AL .75 SH	2 Side Holes	672-035-00		
	AL 1			672-036-00	558-036-00
	AL 1 SH	2 Side Holes	672-037-00		
	AL 1 ST	Short Tip	672-044-00		
	AL 1.5			672-038-00	558-038-00
	AL 2			672-040-00	558-040-00
	AL 2 SH	2 Side Holes	672-041-00		
	AL 2 LBT	Long BRITE TIP™	672-040-0L		
AL 3			672-042-00	558-042-00	

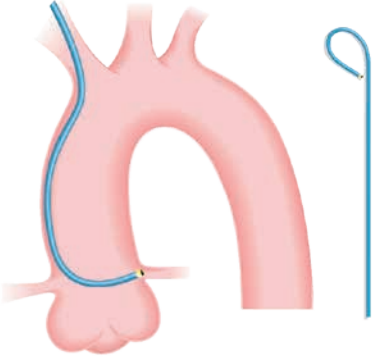

Left Coronary Artery (cont.)	Shape	Description	PRODUCT CODE	
			ADROIT™ 6F	ADROIT™ 5F
 <p>Extra Back-up</p>	XB 2		672-464-00	
	XB 2.5		672-465-00	
	XB 2.75		672-730-00	
	XB 3		672-052-00	558-052-00
	XB 3 SH	2 Side Holes	672-053-00	
	XB 3.25		672-400-00	558-400-00
	XB 3.5		672-054-00	558-054-00
	XB 3.5 (90 cm)		672-054-90	
	XB 3.5 (125 cm)		672-054-25	
	XB 3.5 SH	2 Side Holes	672-055-00	
	XB 3.5	ECOPAK (5 units)	672-054-0E	
	XB 4		672-056-00	558-056-00
	XB 4 SH	2 Side Holes	672-057-00	
	XB 4.5		672-058-00	558-058-00
	XB 4.5 SH	2 Side Holes	672-059-00	

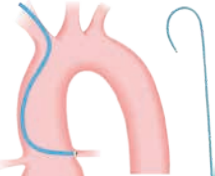
 <p>Extra Back-up LAD</p>	XBLAD 3		672-066-00	
	XBLAD 3 SH	2 Side Holes	672-067-00	
	XBLAD 3.5		672-060-00	558-060-00
	XBLAD 3.5 SH	2 Side Holes	672-061-00	
	XBLAD 4		672-062-00	558-062-00
	XBLAD 4 SH	2 Side Holes	672-063-00	
	XBLAD 4.5		672-064-00	


 <p>Extra Back-up Contralateral</p>	XBC 2.5		672-732-00	
	XBC 2.75		672-734-00	
	XBC 3		672-070-00	558-070-00
	XBC 3 SH	2 Side Holes	672-071-00	
	XBC 3.25		672-736-00	
	XBC 3.5		672-072-00	558-072-00
	XBC 3.5 SH	2 Side Holes	672-073-00	
	XBC 4		672-074-00	558-074-00
	XBC 4 SH	2 Side Holes	672-075-00	
	XBC 4.5		672-076-00	
	XBC 4.5 SH	2 Side Holes	672-077-00	


Right Coronary Artery	Shape	Description	PRODUCT CODE	
			ADROIT™ 6F	ADROIT™ 5F
 <p>Judkins Right</p>	JR 3.5		672-080-00	558-080-00
	JR 3.5 SH	2 Side Holes	672-081-00	
	JR 4		672-082-00	558-082-00
	JR 4 (90 cm)		672-082-90	
	JR 4 (125 cm)		672-082-25	
	JR 4 SH	2 Side Holes	672-083-00	
	JR 4 ST	Short Tip	672-090-00	
	JR 4	ECOPAK (5 units)	672-082-0E	
	JR 5		672-084-00	558-084-00





		Shape	Description	PRODUCT CODE	
				ADROIT™ 6F	ADROIT™ 5F
<b>Right Coronary Artery (cont.)</b>					
Judkins Curved Right		JCR 4		672-098-00	
Amplatz Right		AR 1		672-110-00	558-110-00
		AR 1 SH	2 Side Holes	672-111-00	
		AR 1 LBT	Long BRITE TIP™	672-110-0L	
		AR 2		672-112-00	558-112-00
		AR 2 SH	2 Side Holes	672-113-00	
		AR 2 LBT	Long BRITE TIP™	672-112-0L	
Extra Back-up Right		XBR 1		672-122-00	
		XBR 1 SH	2 Side Holes	672-123-00	
		XBR 2		672-124-00	
Extra Back-up RCA		XB RCA		672-126-00	558-126-00
		XB RCA SH	2 Side Holes	672-127-00	
3 DRC – Williams Technique		3 DRC		672-130-00	558-130-00
		3 DRC SH	2 Side Holes	672-131-00	
NR – Noto Technique		NR 4		672-120-00	
		Shape	Description	PRODUCT CODE	
<b>Dedicated Transradial Shapes</b>				ADROIT™ 6F	ADROIT™ 5F
Jean Fajadet Left		JFL		672-162-00	558-162-00
		JFL ST	Short tip	672-164-00	
Jean Fajadet Right		JFR		672-168-00	
		JFR ST	Short tip	672-170-00	
Radial-brachial Kimny Technique		RB		672-172-00	
Radial Bi-lateral		RBL 4		672-717-00	
		RBL 4.5		672-729-00	
		RBL 5		672-718-00	
Barbeau		BARBEAU		672-174-00	558-174-00
Right Radial		RRAD5			558-320-00
		RRAD6		672-300-00	
Right Radial Extra Back-up		RRAD6 XB		672-301-00	

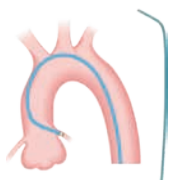
Multipurpose	Shape	Description	PRODUCT CODE	
			ADROIT™ 6F	ADROIT™ 5F
	MPA 1		672-270-00	558-270-00
	MPA 1 (125 cm)		672-270-25	
	MPA 1 SH	2 Side Holes	672-271-00	
	MPA 1 LBT	Long BRITE TIP™	672-270-0L	
	MPB 1		672-272-00	

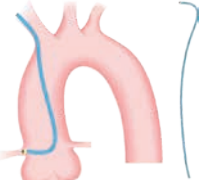	H-STICK		672-278-00	558-278-00
	H-STICK (90 cm)		672-278-90	
	H-STICK SH	2 Side Holes	672-279-00	
	H-STICK LBT	Long BRITE TIP™	672-278-0L	

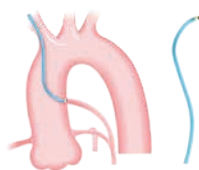Bypass	Shape	Description	PRODUCT CODE	
			ADROIT™ 6F	ADROIT™ 5F
	LCB		672-180-00	558-180-00
	LCB (90 cm)		672-180-90	
	LCB SH	2 Side Holes	672-181-00	

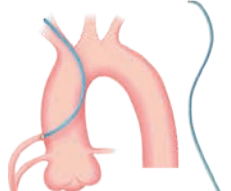	RCB		672-182-00	558-182-00
	RCB (90 cm)		672-182-90	
	RCB SH	2 Side Holes	672-183-00	

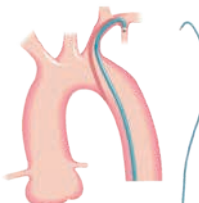	IM		672-190-00	558-190-00
	IM (90 cm)		672-190-90	
	IM SH	2 Side Holes	672-191-00	

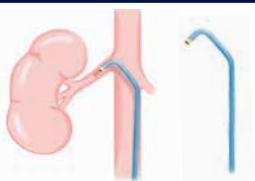Renal Curve	Shape	Description	PRODUCT CODE	
			ADROIT™ 6F	ADROIT™ 5F
	RDC		672-212-55	

## ADROIT™ Guiding Catheter from Cordis

Speak with your Cordis representative about getting more out of your guiding catheter.



For Healthcare Professionals Only.  
 Important information: Prior to use, refer to the instructions for use supplied with this device for indications, contraindications, suggested procedure, warnings and precautions. As part of its continuous product development policy, Cordis reserves the right to change product specifications without prior notification. Please contact your Cordis representative for additional product availability information.  
 CORDIS, Cordis LOGO, ADROIT and VISTA BRITE TIP are trademarks of Cordis and may be registered in the US and/or in other countries.  
 All other marks are the property of their respective owners. © 2023 Cordis. All Rights Reserved. 100528428-5 11/2023

