


# ILLUMINATING THE PATH TO PRECISION PLACEMENT



**BRITE TIP™**  
INTERVENTIONAL SHEATH INTRODUCER

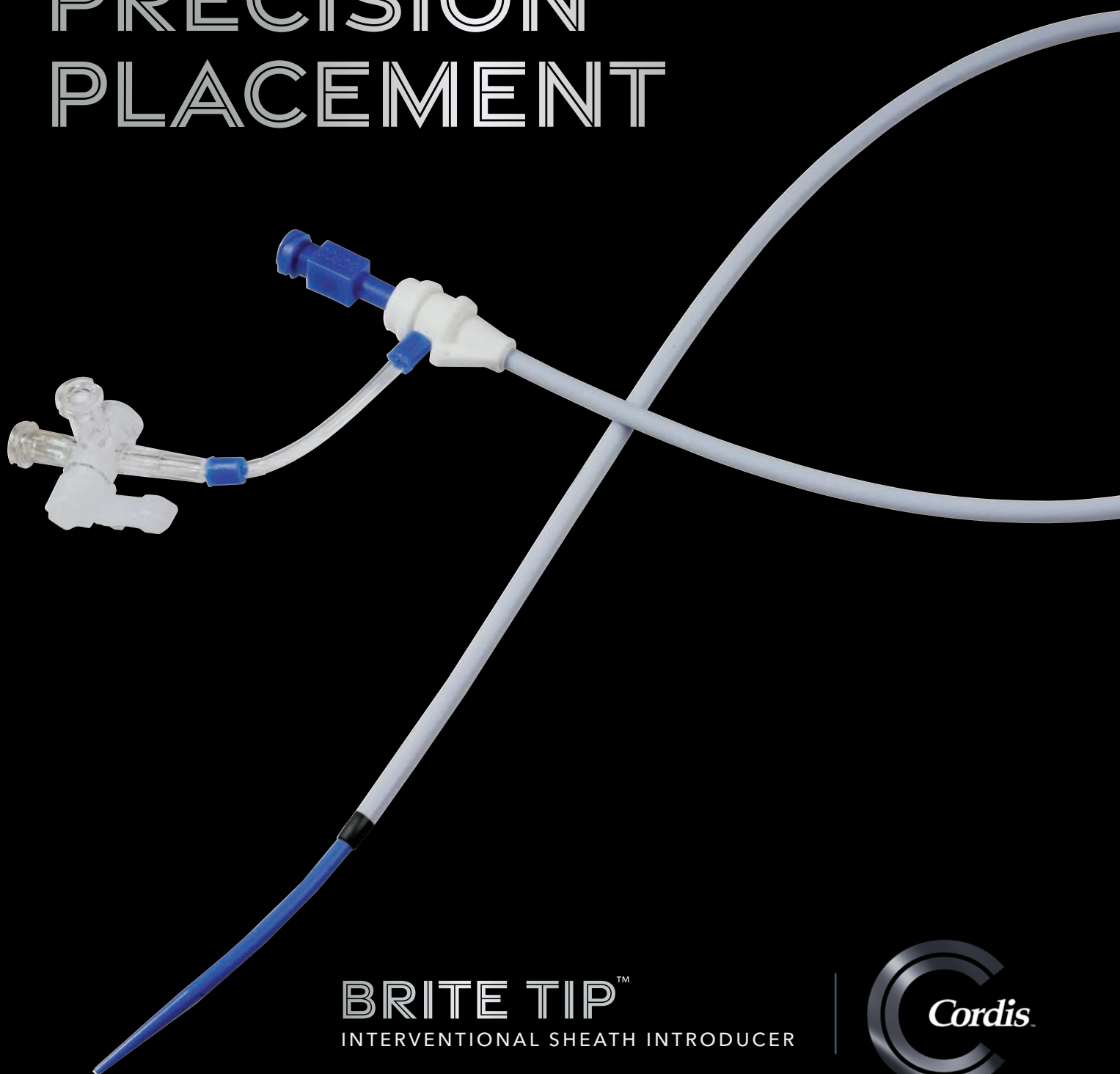

# KINK RESISTANT CANNULA

## RESILIENT DUAL LAYER CANNULA


- Flexibility and support to access tortuous vessels

## REINFORCED CANNULA WALL

- Kink resistance

## SILICONE COATED CANNULA AND VESSEL DILATOR

- Lubricious, anti-thrombogenic vessel entry

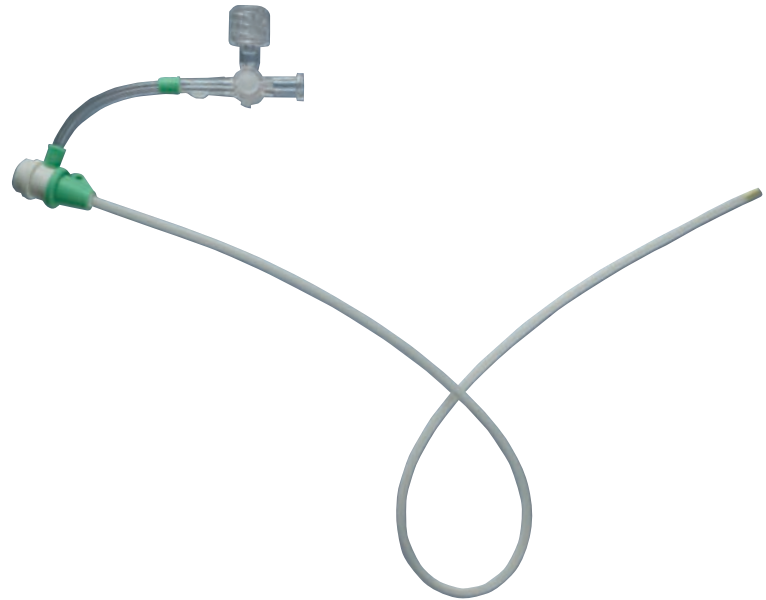

# SLIX VALVE

## SIX-CUT SPIRAL DESIGN FOR HEMOSTATIC VALVE

- Free catheter movement without blood reflux or air aspiration




# **BRITE TIP™** Sheath Introducers

Precise, atraumatic placement

## **MANICURED TIP AND DILATOR DESIGN**

### **TAPERED TRANSITIONS ON SHEATH AND DILATOR FOR:**

- Atraumatic vessel entry
- Improved patient comfort

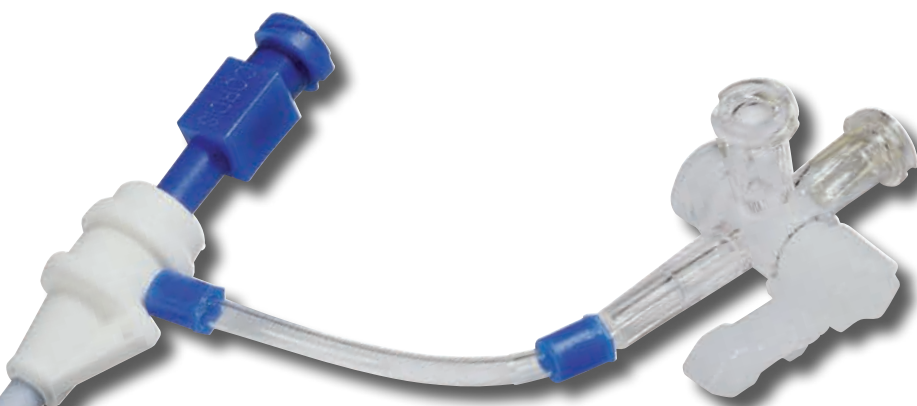

## **RADIOPAQUE BRITE TIP™ SHEATH**

### **BLACK TUNGSTEN-FILLED DISTAL TIP**

- Excellent visualization facilitates exact location of sheath's distal tip for accurate interventional positioning



# BRITE TIP™

## Sheath Introducers

### ORDERING INFORMATION

Longer lengths available for a wide variety of interventional procedures



### CATALOG NUMBERS

Cannula Length			4.5cm	5.5cm	11cm	23cm	35cm	45cm	55cm	70cm	90cm
French Size	Color Code	Dilator Tip (mm)	With Mini-Wire	With Mini-Wire	With Mini-Wire	Without Mini-Wire	Without Mini-Wire	Without Mini-Wire	Without Mini-Wire	Without Mini-Wire	Without Mini-Wire
4	■	25	SRA0116		401411M	401423M	SRA0207	SRAM17661			
5	■	25	SRA0117	SRA0035	401511M	401523M	401535M	401545M	401555M		
6	■	35	SRA0115	401605M	401611M	401623M	401635M	401645M	401655M	SRA0217	401690M
7	■	35		401705M	401711M	401723M	401735M	401745M	401755M		401790M
8	■	35		401805M	401811M	401823M	401835M	401845M	401855M		
9	■	35			401911M	401923M	401935M				
10	■	45			401011M	401023M					
11	■	45			401111M	401123M					

\*Each catalog number must be ordered by the box

### QUANTITY

- 5 units per package for 4.5cm, 5.5cm, 11cm and 23cm lengths
- 2 units per package for 35cm or longer lengths

For Healthcare Professional Only. Important information: Prior to use, refer to the instructions for use supplied with this device for indications, contraindications, suggested procedure, warnings and precautions. As part of its continuous product development policy, Cordis reserves the right to change product specifications without prior notification.

Please contact your Cordis representative for additional product availability information.

CORDIS, Cordis LOGO and BRITE TIP are trademarks of Cordis and may be registered in the US and/or in other countries.  
© 2026 Cordis. All Rights Reserved. 100682856

