


TEMPO AQUA[™]

HYDROPHILIC-COATED DIAGNOSTIC CATHETER

Navigating most challenging pathways with lubricity and control

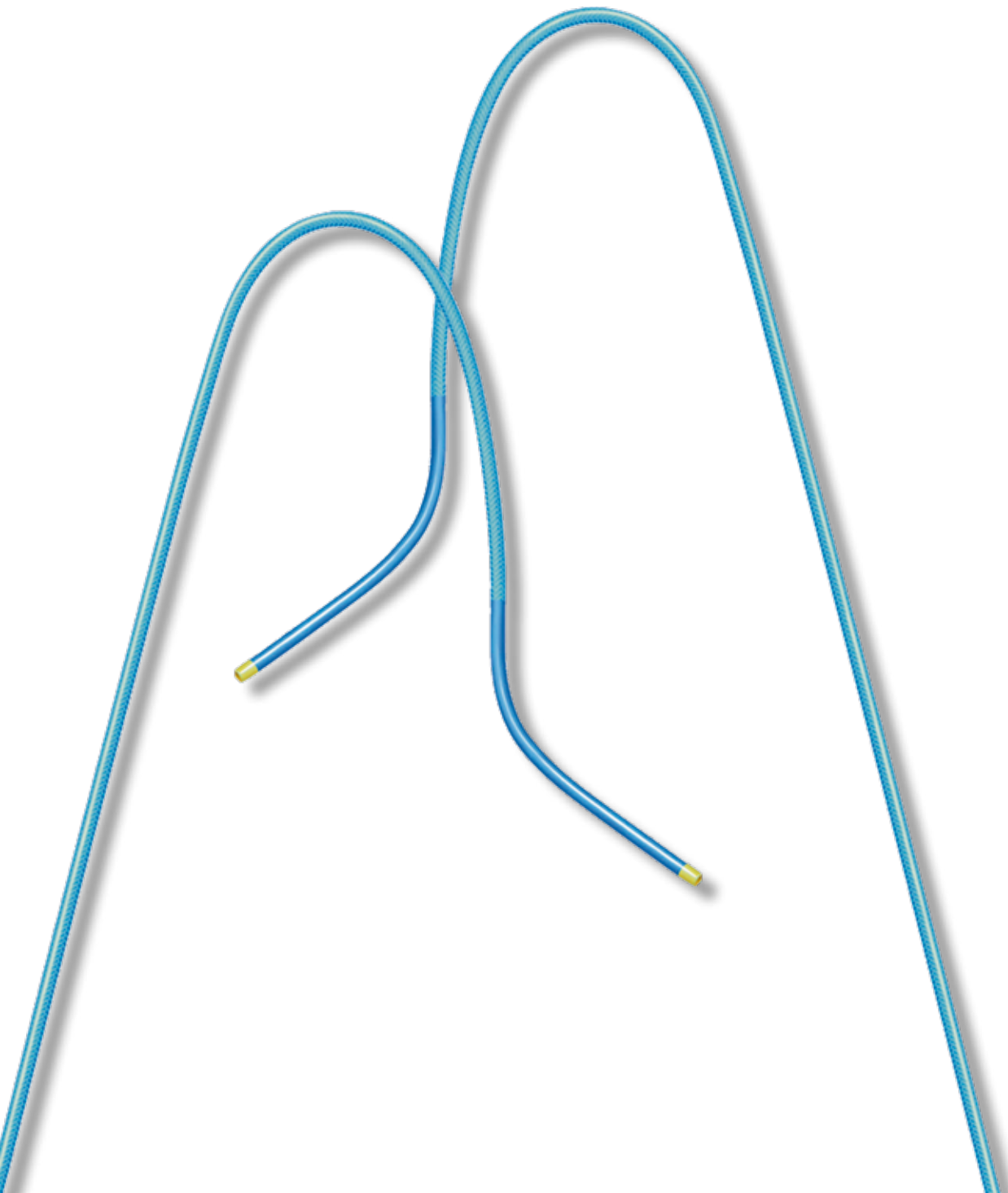

The **TEMPO AQUA™** Catheter makes endovascular procedures quick and easy by combining the Cordis TEMPO™ Angiographic Catheter characteristics with a durable hydrophilic coating for easy advancement into the vasculature and subintimal space.

CROSSABILITY

- Hydrophilic coating to **navigate the most challenging cases** with lubricity and control
- Hydrophilic coated distal shaft enables **smooth passage through the vessels**
- Catheter tip demonstrates superb shape retention, greater positional stability and is made from a special blend of soft **atraumatic** material

CONTROL

- Distal half of the catheter is coated for high lubricity, **proximal part is uncoated** for secure handling
- High-strength nylon construction and braided shaft design give strength and flexibility for **superb torque response** and maximum **tactile control**
- Thin-wall technology guarantees a large inner diameter and provides **high flow rates at low injection force**

SELECTIVE CATHETER ORDERING INFORMATION

TEMPO AQUA™ Catheter - 4 French

Catalog Number	Description	Shape	Side Holes	Length (cm)
452-442V0	Cobra I	C1	0	65
452-443V0	Cobra II	C2	0	65
SRD6913	Straight	STR	0	65
SRD6912	Vertebral 135°	VER	0	65
452-443S0	Cobra II	C2	0	80
SRD6908	Straight	STR	0	90
452-413H0	Berenstein	BER	0	100
452-435H0	Headhunter I	H1	0	100
452-407H0	Multipurpose A-1, normal curve	MPA-1	0	100
452-430H0	Sidewinder I	SIM1	0	100
452-431H0	Sidewinder II	SIM2	0	100
SRD6821	Straight	STR	0	100
452-414H0	Vertebral 135°	VER	0	100
SRD6107	Headhunter I	H1	0	125
SRD6667	Straight	STR	0	125
SRD6108	Vertebral 135°	VER	0	125

TEMPO AQUA™ Catheter - 5 French

Catalog Number	Description	Shape	Side Holes	Length (cm)
452-542V0	Cobra I	C1	0	65
452-543V0	Cobra II	C2	0	65
SRD6658	PIER	PIER	0	65
SRD6906	Vertebral 135°	VER	0	65
SRD6677	Cobra II	C2	0	80
SRD6690	Modified Headhunter I	MOD H1	0	90
452-513H0	Berenstein	BER	0	100
452-535H0	Headhunter I	H1	0	100
452-507H0	Multipurpose A-1, normal curve	MPA-1	0	100
452-530H0	Sidewinder I	SIM1	0	100
452-531H0	Sidewinder II	SIM2	0	100
452-514H0	Vertebral 135°	VER	0	100
SRD6905	Straight	STR	0	100
SRD6721	Headhunter I	H1	0	125
SRD6109	Sidewinder II	SIM2	0	125
SRD6669	Straight	STR	0	125

For Healthcare Professionals Only.

Important information: Prior to use, refer to the instructions for use supplied with this device for indications, contraindications, side effects, suggested procedure, warnings and precautions. As part of its continuous product development policy, Cordis reserves the right to change product specifications without prior notification. Please contact your Cordis representative for additional product availability information.

CORDIS, Cordis LOGO, TEMPO are trademarks of Cordis and may be registered in the US and/or in other countries. © 2021 Cordis. All Rights Reserved. 100589211 08/2021